



Distrattori
Cranio - maxillo - facciali
Cranio - maxillofacial distractors

Descrizione

Apparecchio intraorale per la distrazione controllata e misurabile del ramo/corpo della mandibola:

- luxtaosseo con asta di attivazione intraorale;
- Di facile applicazione sia a destra che a sinistra modellando e tagliando le placche;
- Disponibile in tre versioni con ampiezza di distrazione fino 15 mm, 20 mm e 30 mm.

Vantaggi

- Il distrattore è molto sottile, inferiore a 4 mm, e, quindi, poco avvertibile dal paziente.
- Di facile posizionamento, perché l'ampiezza della placca consente di mettere le viti dove è più favorevole.
- Applicabile a destra e a sinistra. Riduce lo stock di magazzino.
- Facilità di adattamento alla conformazione ossea. Si adatta a tutte le mandibole.
- Possibilità di togliere, ritagliandola, la parte di placca che non serve.
- A distrazione ultimata si rimuove la bussolina intraorale e, in attesa del consolidamento del rigenerato osseo, il distrattore rimane tutto sotto mucosa.
- Da impiantare con viti di Ø 2,0 Cizeta Surgical e relativo strumentario.

Description



Intraoral device for the controlled and measurable distraction of the mandibular ramus/body:

- *luxtaosseous with intraoral activation bar;*
- *Easily applicable on both left and right side on mandible by modeling and cutting the plates;*
- *Available in three patterns with distraction range up to 15 mm, 20 mm and 30 mm.*

Benefits

- *Low profile distractor, less than 4 mm thickness, hardly noticeable by the patient.*
- *Easily positioned, since the plate's width allows screws to be placed in the most favourable position.*
- *It can be applied on both left and right side of mandible. It minimizes warehouse stock.*
- *Easily adapted to bone conformation. It adapts to all mandible conformation.*
- *The part of plate in excess can be removed by cutting it.*
- *At the end of distraction period the intraoral thimble is removed and, waiting for the consolidation of the regenerate bone, the distractor remains entirely under the mucous membrane.*
- *To be implanted with Ø 2,0 Cizeta Surgical screws and related instruments.*

01 - CAPE Distrattore del ramo/corpo mandibolare

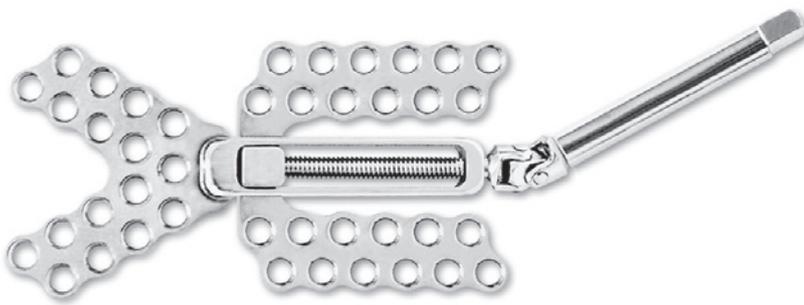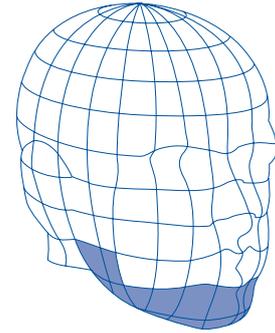

01 - CAPE Mandibular ramus/body distractor

Cranio - maxillofacial distractors